



Ring Placement

1

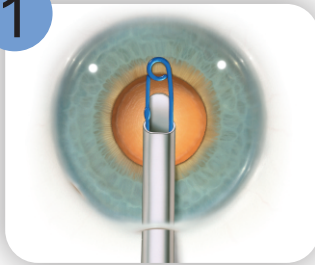


With viscoelastic in the Anterior Chamber inject the Ring approximately 2mm forward, then engage the distal iris margin.

2

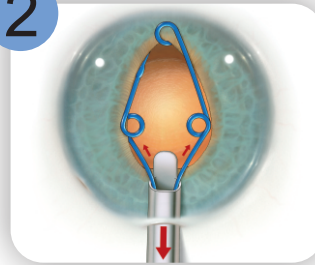


Inject while simultaneously retracting the injector. "Toeing" the injector down eases placement of the lateral scrolls.

3

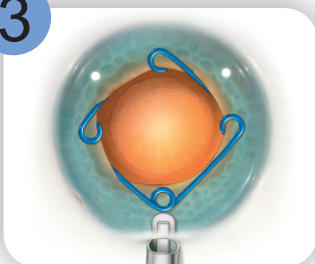

Lay proximal scroll on anterior iris.


4

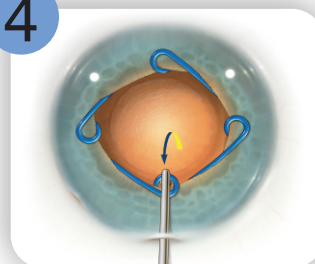

Use a Osher/Malyugin Ring Manipulator to place the proximal scroll.


Ring Removal

1

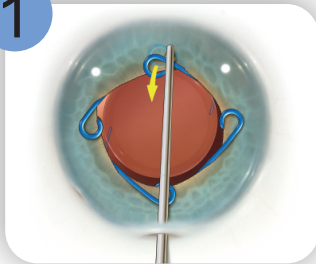


Use a Osher/Malyugin Ring Manipulator to disengage the distal scroll.

2

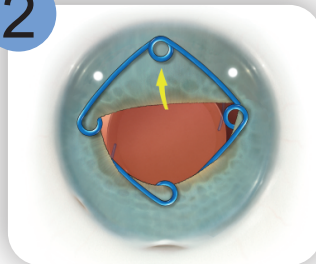


Place the distal scroll onto anterior iris.

3

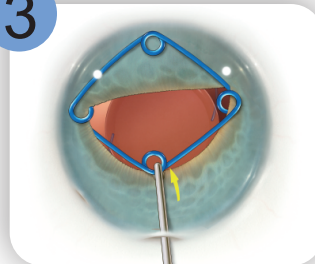


Move the proximal scroll toward the center of the pupil, and a bit anterior to disengage.

4

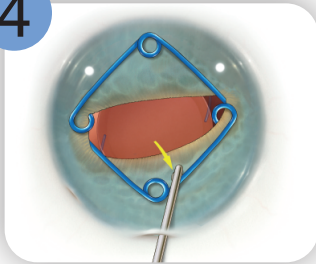


Place the proximal scroll onto anterior iris.

5

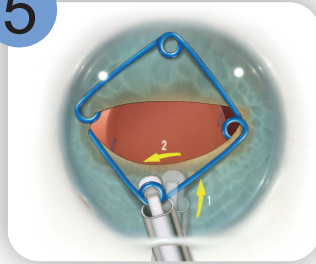

Position the injector's platform on the right under the ring, move under the proximal scroll, then move the injector button forward.

6

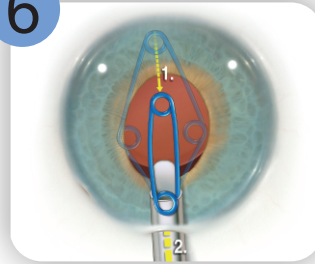

Move the injector button back to remove the ring. Stop when the left scroll rides above the cannula, then remove the "assembly".